

Field Service Procedure

Double Vapor Parking Accessory Anesthesia



Revision -
12/19/03
SP00301

[RETURN TO SERVICE PROCEDURE TABLE OF CONTENTS](#)
[RETURN TO CD-ROM TABLE OF CONTENTS](#)

Service Procedure

Double Vapor Parking Accessory

Installation of Double Vapor Parking Accessory

1 PURPOSE

The purpose of this document is to provide instruction for the installation of the Vapor Parking Accessory onto the anesthesia machine trolley.

The following parts are contained with the kit:

1. Double Vapor Parking Mount (P/N's listed in [Fig. 1](#))
2. Stop Block (P/N 4115872)

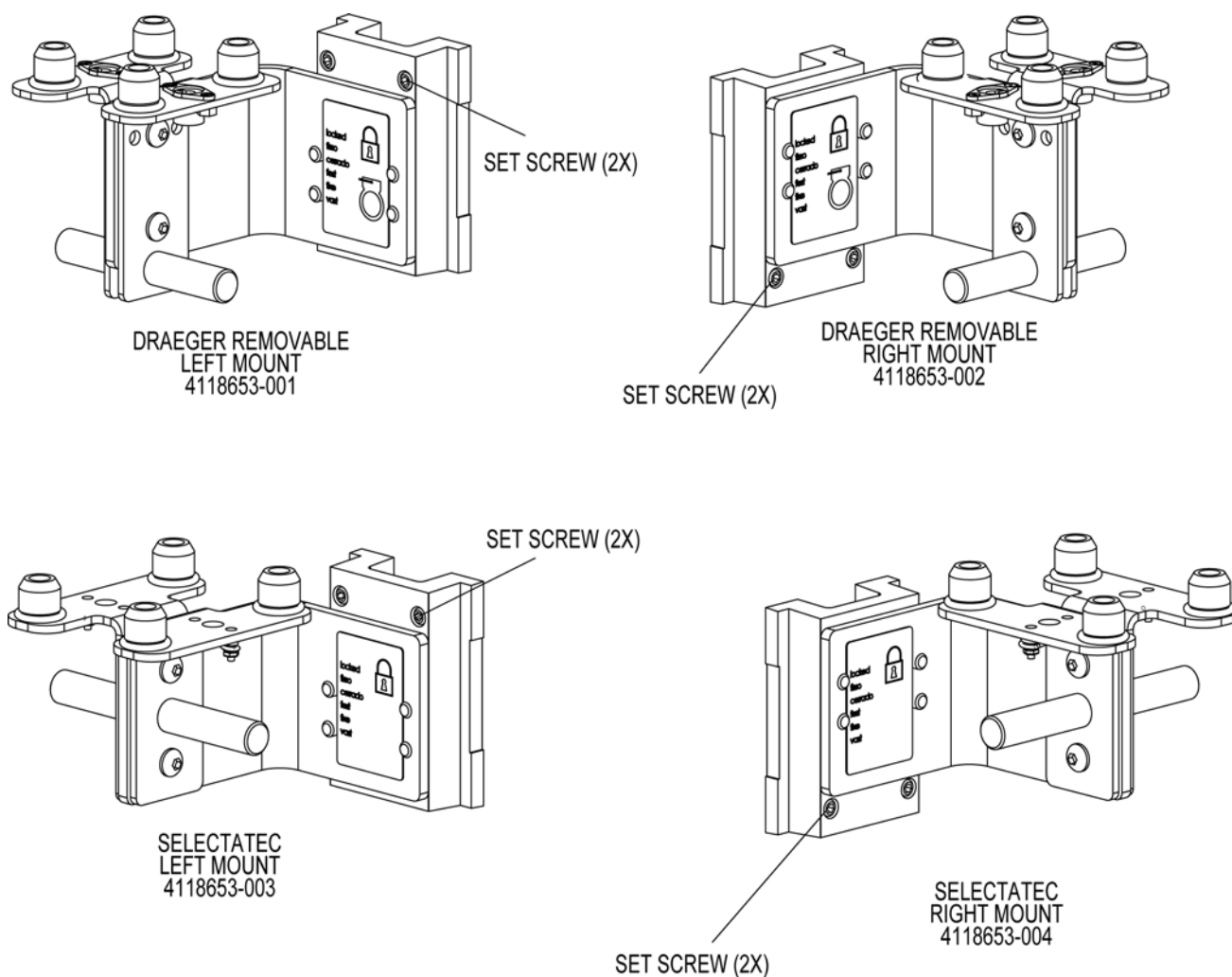


Fig. 1 Double Vapor Parking Mount Assemblies

Double Vapor Parking Accessory

Service Procedure

2 INSTALLATION PROCEDURE

1. The Double Vapor Parking Mount will be installed on the trolley on the opposite side from the breathing system.
2. Install the Stop Block into the GCX rail.



WARNING:

Failure to install the Double Vapor Parking Mount within the specified range can lead to system instability and result in equipment damage or personal injury.

3. Install the Double Vapor Parking Mount into the GCX channel. Position the mount at the height indicated in Fig. 2 or Fig. 3. Tighten the two 3 mm set screws against the channel. Refer to Fig. 1 for set screw locations.
4. Position the Stop Block against the bottom of the Double Vapor Parking Mount and tighten the Stop Block Screw against the channel.

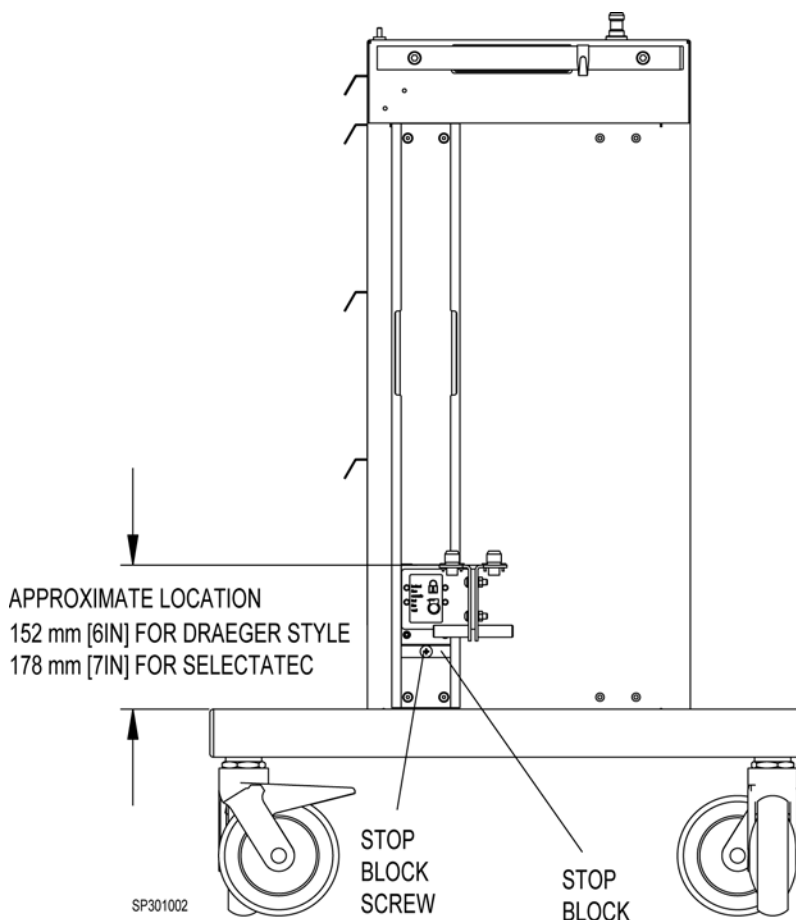


Fig. 2 Installation on Right Side of Machine

Service Procedure

Double Vapor Parking Accessory

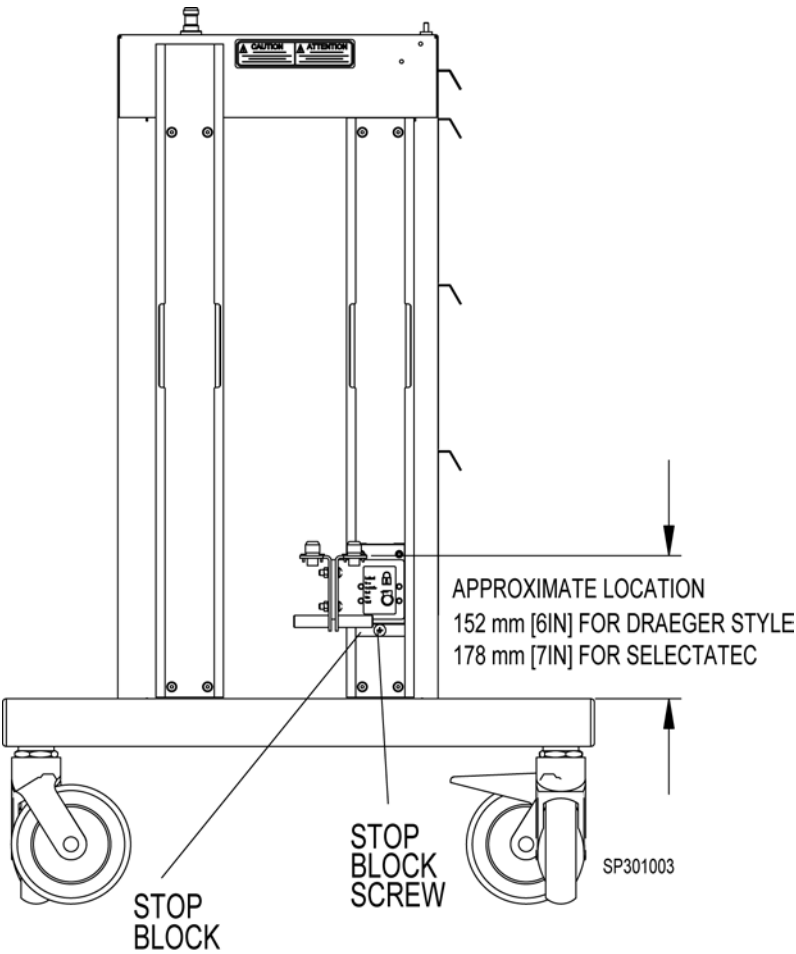


Fig. 3 Installation on Left Side of Machine

Double Vapor Parking Accessory

Service Procedure

In the United States:
Draeger Medical, Inc.
3136 Quarry Road
Telford, PA 18969
USA

Tel: (215) 721-5402
(800) 543-5047
Fax: (215) 721-5785

Double Vapor Parking Accessory en
SP00301

Subject to change without notice.
Will not be replaced in the event of modifications.

In all other countries:
Dräger Medical AG & Co. KGaA
Moislinger Allee 53 – 55
D-23542 Lübeck
Germany

Tel: (++49) (0) 1805-3723437
Fax: (++49) 451/882 - 3779